

## Your image, when and where you need it

POLYMOBIL Plus, POLYMOBIL III

Answers for life.

**SIEMENS**



## POLYMOBIL Plus, POLYMOBIL III

# Simplicity and reliability in mobile X-ray

Over the years, and in thousands of facilities worldwide, POLYMOBIL® has established its strong reputation. The POLYMOBIL family of mobile generators lives up to its name in offering mobile versatility for a multitude of applications. The Siemens POLYMOBIL III and POLYMOBIL Plus offer you the choice between a “basic” version and a “basic with extras”. POLYMOBIL Plus boasts a 20 kW generator with 300 mAs – more power for higher image quality.

POLYMOBIL makes it easy to bring the X-ray system to the patient. Whether for patient wards, intensive care, operating room, sports medicine practice, occupational medical facilities or for veterinary use, POLYMOBIL is designed to give you high value for the investment.





SIEMENS

POLYMOBIL Plus

## POLYMOBIL Plus, POLYMOBIL III

# Mobile X-ray system directly at the patient's bedside

### Easy to move

Transporting POLYMOBIL from department to department can include the necessity to move it in and out of elevators, over thresholds, and around corners. The platform's short wheelbase with large polyurethane tires helps to overcome rolling resistance, ensuring that the system glides easily.

The tube support arm is counterbalanced, allowing great freedom of movement and precise positioning. The adjustable collimator and light marker ensure rapid and correct positioning in the desired projection. Because the tube head is rotatable in all planes, both horizontal and vertical exposures can be made – and the X-ray tube angulation is conveniently shown on the scales on both sides of the tube head.


### Easy to use

POLYMOBIL Plus and POLYMOBIL III are as easy to use as they are to move. The touch-screen keys and digital display of kV and mAs values on the ergonomically designed operating console enable exposure parameters to be set easily and quickly.

The high-frequency, high-output X-ray generator with multipulse voltage provides high dose yield, keeping exposure times short and minimizing motion artifacts. This enables good reproducibility and results in consistent image quality.

### POLYMOBIL easily meets your expectations.

Just plug POLYMOBIL into any standard wall socket and you have full power – always.

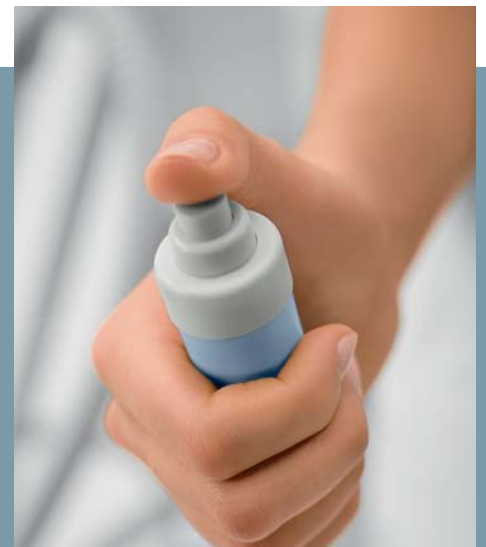

These features put POLYMOBIL right at the top of the list for basic mobile X-ray:

#### Easy to choose

When system reliability, simple features for streamlined operation, and value-for-money quality are your priorities for mobile X-ray, POLYMOBIL from Siemens is your choice.

With over 100 years of experience in medical imaging, Siemens is known worldwide for its top-quality clinical solutions, from basic mobile radiographic technology to the broadest spectrum of complex dedicated diagnostic and therapeutic systems. With a worldwide distribution network in over 140 countries, Siemens offers equipment sales, financing, and service.

Experience the versatility and simplicity of POLYMOBIL in your facility.

- Transport the light-weight and compact POLYMOBIL easily within the facility
- Short wheelbase with large polyurethane tires helps to glide easily
- Use the simple digital display for setting the two exposure parameters – kV and mAs
- Set the exposure field precisely with the easy-to-operate collimator
- Move the counterbalanced tube arm and rotatable tube head for precise adjustment – effortlessly



POLYMOBIL Plus and POLYMOBIL III are not commercially available in the United States (U.S.A.).

The information in this document contains general descriptions of the technical options available and may not always apply in individual cases.

The required features should therefore be specified in each individual case at completion of contract. Siemens reserves the right to modify the design and specifications contained herein without prior notice. Please contact your local Siemens sales representative for the most current information.

Original images always lose a certain amount of detail when reproduced.

In the interest of complying with legal requirements concerning the environmental compatibility of our products (protection of natural resources, waste conservation), we recycle certain components. Using the same extensive quality assurance measures as for new components, we guarantee the quality of these recycled components.

Please find fitting accessories:  
[www.siemens.com/medical-accessories](http://www.siemens.com/medical-accessories)

## Local Contact Information

### In Japan

Siemens-Asahi  
Medical Technologies Ltd.  
Takanawa Park Tower  
3-20-14, Higashi-gotanda  
Shinagawa-ku  
Tokyo 141-8644  
Phone: +81 3 5423 8510

### In Asia

Siemens Pte Ltd  
The Siemens Center  
60 MacPherson Road  
Singapore 348615  
Phone: +65 6490 8182

### In Germany

Siemens AG  
Healthcare Sector  
Karlheinz-Kaske-Str. 2  
91052 Erlangen  
Phone: +49 9131 84-0

## Global Business Unit

Siemens AG  
Medical Solutions  
Special Systems  
Henkestr. 127  
DE-91052 Erlangen  
Germany  
Phone: +49 9131 84-0

## Global Siemens Headquarters

Siemens AG  
Wittelsbacherplatz 2  
80333 Muenchen  
Germany

## Global Siemens Healthcare Headquarters

Siemens AG  
Healthcare Sector  
Henkestr. 127  
91052 Erlangen  
Germany  
Phone: +49 9131 84-0  
[www.siemens.com/healthcare](http://www.siemens.com/healthcare)

## Legal Manufacturer

Siemens AG  
Wittelsbacherplatz 2  
DE-80333 Muenchen  
Germany

[www.siemens.com/healthcare](http://www.siemens.com/healthcare)